



Field Service Procedure

Part Number: SP00275

Rev: A

Date: 18 February 2004

© 2004 Draeger Medical, Inc.

Domestic Suction Mount on Fabius Tiro™

[RETURN TO SERVICE PROCEDURE TABLE OF CONTENTS](#)
[RETURN TO CD-ROM TABLE OF CONTENTS](#)

GCX Mounting Assembly Installation Guide

Purpose: The purpose of this document is to provide instruction for the mounting of Suction Assembly Mounting Kit to Fabius Tiro Anesthesia Machines.

PARTS REFERENCE:

The following parts and hardware are included with this installation kit (see illustrations for parts, hardware not shown):

Item #	Description	Qty
1	Rail Assembly	1
2	Regulator Mounting Bracket	1
3	Canister Mounting Bracket	1
4	Canister Hoop	1

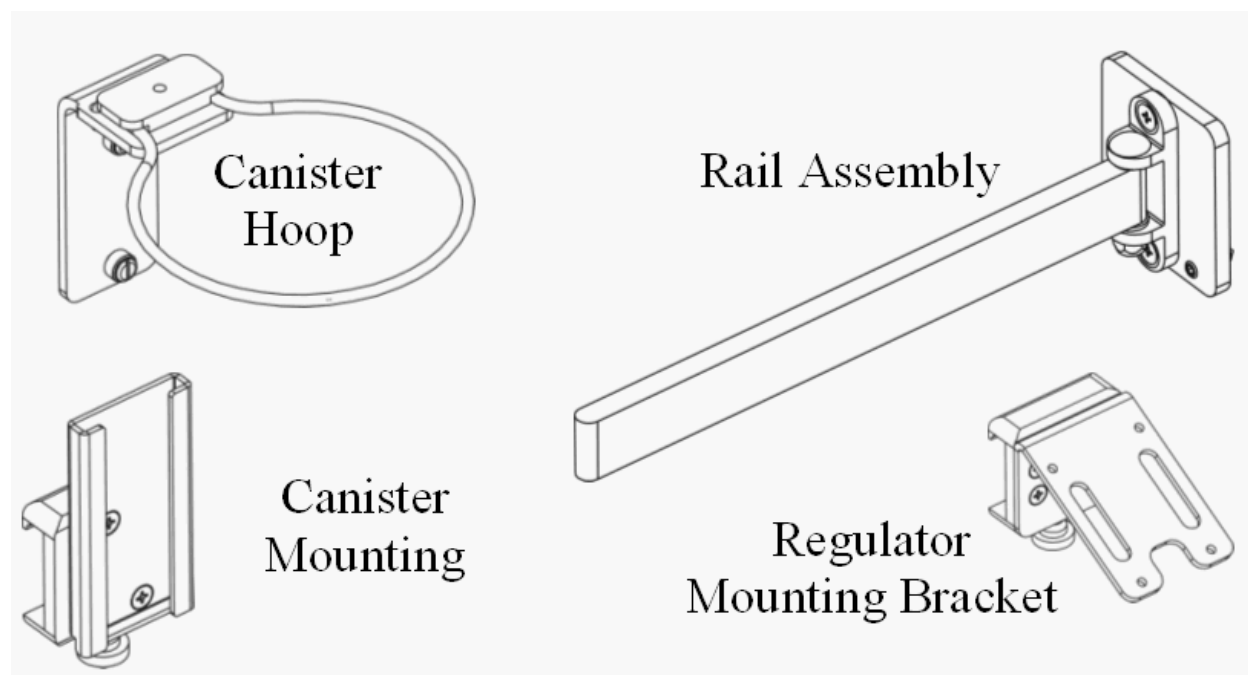


Figure 1. GCX Suction Assembly

Tools Required for This Procedure: Phillips screwdriver, 3mm hex key (provided).

SP00275
9/18/03

INSTALLATION PROCEDURE:

Mounting the Rail Assembly

1. Insert the Slide Rail Assembly (w/ screws) into cutout in GCX channel and slide the Rail Assembly to the desired height.
2. Using the 3mm hex key (provided), tighten the two (2) clamping screws to secure the Assembly in the channel.

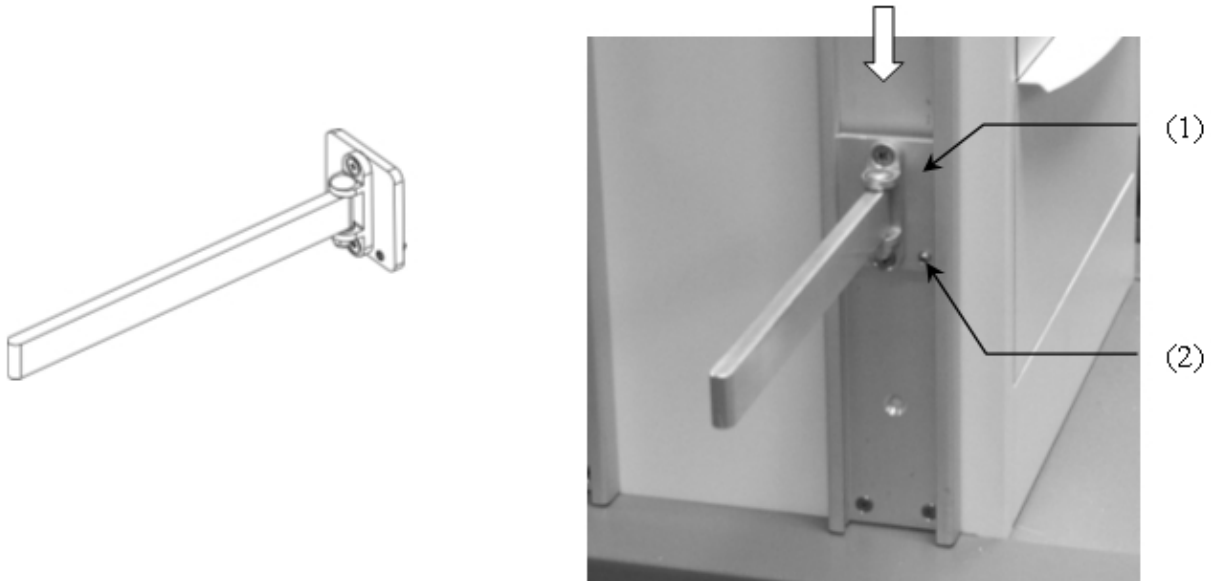


Figure 2. Rail Assembly

Mounting the Regulator

1. Unscrew four (4) existing screws (ITEM 1) from the bottom of the regulator (shown below). Retain screws.
2. Attach the Regulator Mounting Bracket (ITEM 2) to the regulator using the four (4) screws removed in Step 1.

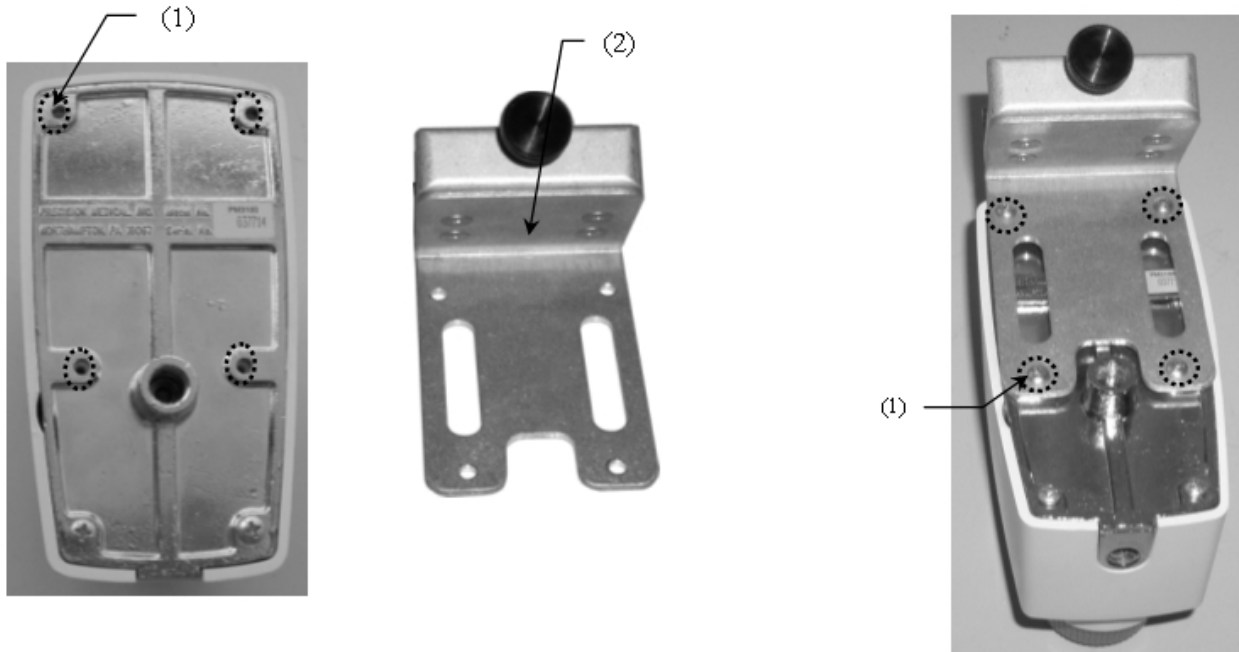


Figure 3. Mounting Regulator Bracket

3. Hook the Regulator Mounting Bracket over the rail. Tighten the clamping knob (bottom of Mounting Bracket) to secure the regulator to the rail.



Figure 4. Rail Assembly

SP00275
9/18/03

Mounting the Hose Holder

1. Hook the Hose Holder over the 10mm rail. Push down to secure the Holder to the rail.

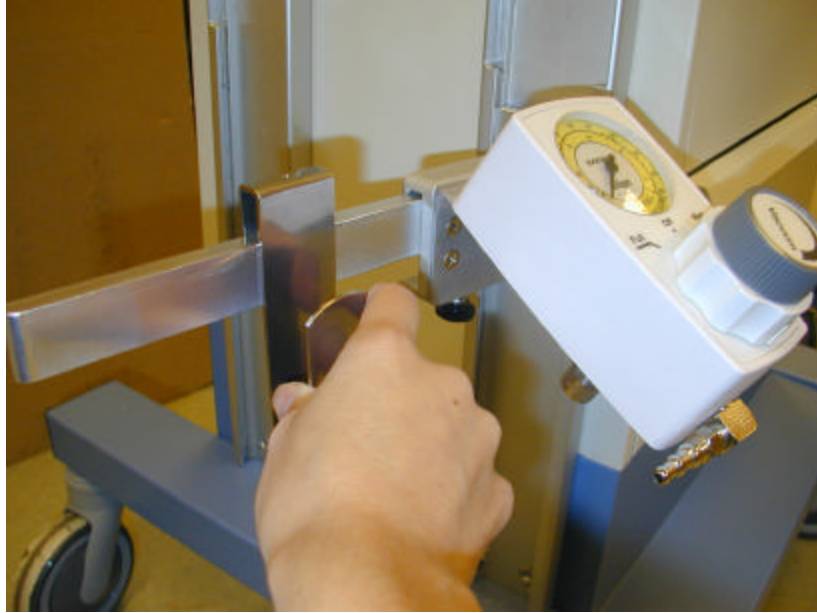
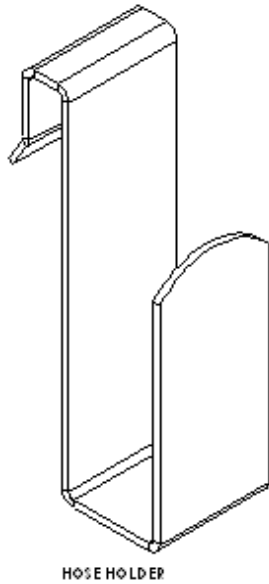


Figure 5. Mounting Hose Holder

Mounting the Suction Canister Bracket

1. Hook the Canister Mounting Bracket over the 10mm rail. Tighten the clamping knob (bottom of Mounting Bracket) to secure the Bracket to the rail.

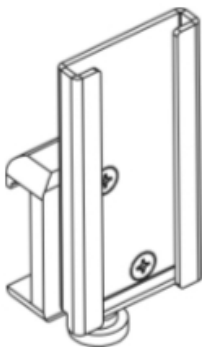


Figure 6. Mounting Suction Canister Bracket

2. Mount the suction canister (with mounting detail) in the channel (not shown).
Installation Note: If the suction canister does not provide a mounting detail, see *"Mounting the Canister Hoop."*

Mounting the Canister Hoop

1. Slide the Canister Hoop into the channel and tighten the two Nylon screws.
2. Insert the suction canister in the Canister Hoop.

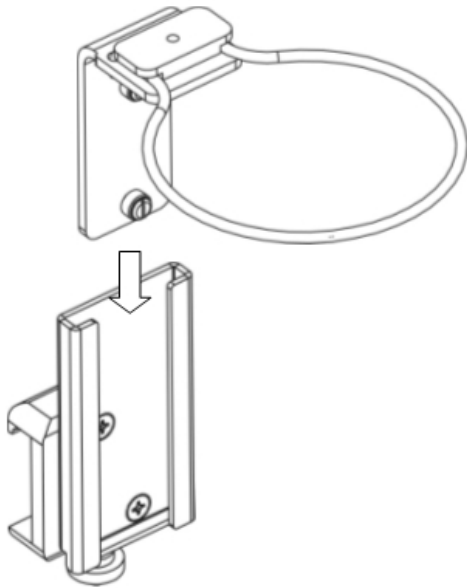


Figure 7. Mounting Canister Hoop

[RETURN TO SERVICE PROCEDURE TABLE OF CONTENTS](#)
[RETURN TO CD-ROM TABLE OF CONTENTS](#)

[RETURN TO SERVICE PROCEDURE TABLE OF CONTENTS](#)
[RETURN TO CD-ROM TABLE OF CONTENTS](#)



DrägerService is a division of
Draeger Medical, Inc.
3122 Commerce Drive
Telford, PA 18969
Tel: (215) 721-5402
(800) 543-5047
Fax: (215) 721-5784
Web: www.draegermedical.com
Printed in the U.S.A.